

Technical

AM SCIENTIFICS

Bulletin



Dear Surgeon,

Cymetra, Micronized AlloDerm Tissue, is a non-surgical product that restores the defects in soft tissue that cause wrinkles, creases, scars and other signs of aging. **Cymetra** is available upon order.

Cymetra begins as natural tissue that contains no human cells, yet retains collagen and elastin patterns that structurally mimic normal, healthy skin. Once implanted by needle into treatment area, **Cymetra** repopulates with patients own cells - actually reconstructing his/her own skin.

Cymetra looks, feels, and ages exactly like the patient's own natural tissue.

Simply put, **Cymetra** helps restore what nature has taken away. The new tissue supports the patient's skin the way his/her original soft tissue did, eliminating creases, sagging, or depressed scars.

It is important to remember that because **Cymetra** is true tissue replacement, patients vary in how their bodies react to treatment. It is possible that to achieve the desired goal, more than one treatment may be necessary to achieve adequate tissue accumulation. Therefore, your treatment plan should begin with an in-depth consultation to communicate and evaluate the desires and expectations of your patient and how you would proceed.

Cymetra is commonly applied to:

- **DEPRESSED SCARS** - Reduces the appearance of depressed scars from surgery, injury, or disease, including acne scars.
- **LIP RESTORATION** - Restores lips to their original potential.
- **NASOLABIAL FOLDS** - Diminishes the creases that run from the corner of the mouth to the nose, giving face a smoother overall look.
- **FACIAL CREASES** - Smoothes deep wrinkles and facial creases, including marionette and laugh lines.

Safety First

Cymetra is implanted through a simple procedure that can be performed right in your physician's office. Each application includes a pain-relief medication to ensure maximum comfort.

Cymetra is a form of **AlloDerm** Tissue - a natural skin product with no reported viral transmission in 12 years of use and more than 450,000 grafts. Because human body will not recognise **Cymetra** as a foreign substance, there's minimal risk of tissue rejection or allergic reaction. In scientific studies, no cases of allergic reaction to **Cymetra** have been reported.

AM SCIENTIFICS LTD is the exclusive Distributor of **AlloDerm**® in Switzerland. We assure you of our very best services always,

Yours truly,


A.H. MOGHBEL
Director



Cymetra™
Micronized AlloDerm™ Tissue


Recommended Treatment for Atrophic Lips

Vermillion Border (white roll)

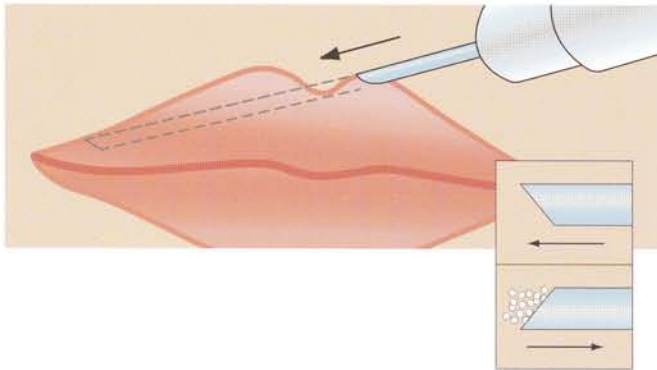


22-23 gauge, 1.5 inch needle

Insert bevel down needle
at peak of cupid's bow,
threading across to the
opposite commissure

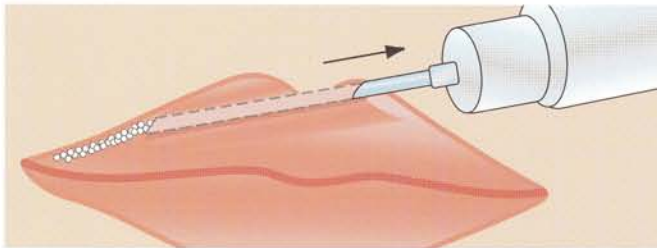

Rotate bevel up and inject
while withdrawing...

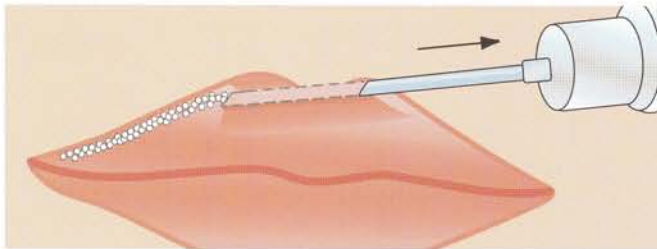


Tent the skin with needle to
ensure superficial placement

Repeat on other side, criss-
crossing the middle third of lip
for double amount of Cymetra

0.5-0.6 units per side

2-3 treatments over time is
recommended for a progressive
buildup of tissue

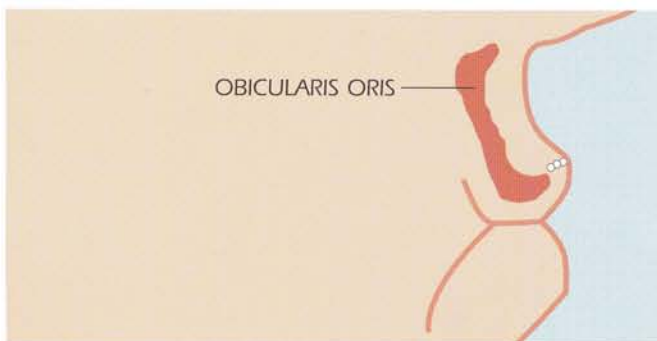

PLANE OF PLACEMENT


Superficial along the white roll



Cymetra™
Micronized AlloDerm® Tissue

Recommended Treatment for Atrophic Lips


Vermillion Area

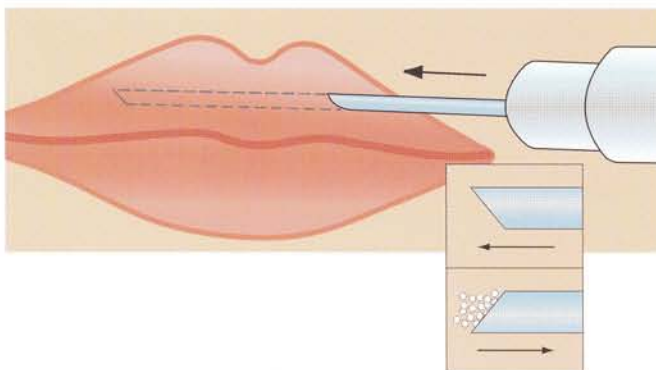


23 gauge, 1.5 inch needle

Insert bevel down needle deep, just anterior to fascia of obicularis muscle one third across the lip

Thread across to the opposite side and **inject with bevel up** while withdrawing

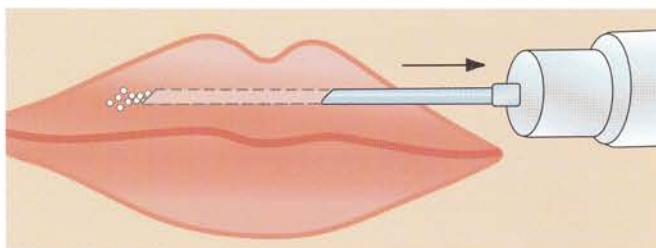

Repeat from the other side, forming a criss-cross pattern in the middle third to deliver double Cymetra in mid-lip

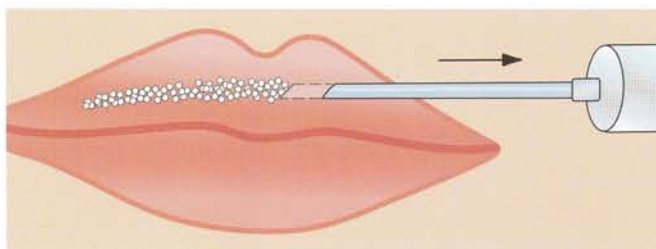


Mold to symmetry, roll anteriorly

No patient massage

No movement of lips for 6 hours or more

1-2 units per lip

2-3 treatments over time is recommended for a progressive buildup of tissue

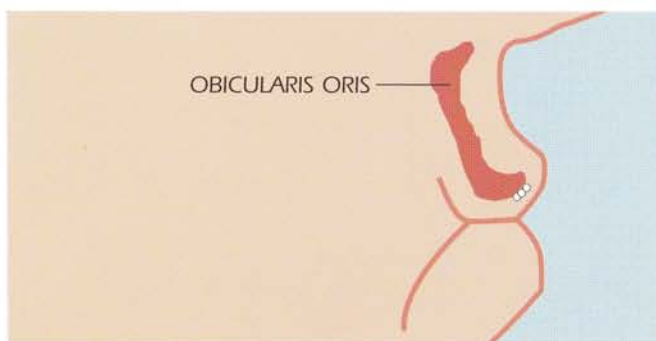

PLANE OF PLACEMENT
Inject along the dermal subcutaneous junction



Established in 1984, AM SCIENTIFICS LTD is active in marketing and distribution of high quality wound care, medical skin care devices and implants destined for use in hospitals, clinics and medical practices.

AM SCIENTIFICS is Exclusive Swiss Distributor for **AlloDerm**®, **Cymetra**™, **JUC**™ and other plastic, burn surgery & dermatology products and is active in most European Union (EU) countries.



AM SCIENTIFICS LIMITED

P.O. Box 507
CH -1212 Grand-Lancy 1
Geneva – Switzerland
T:+41 (0)223011794, F:+41 (0)223011795
M: +41 (0)763372100
www.am-scientifics.com
marketing@am-scientifics.com